


hCG Rapid Pregnancy Test

QUICK REFERENCE GUIDE



1. Allow the Pregnancy Test device and urine sample to reach room temperature 15-30°C (59-86°C) before opening the foil pouch.
2. Remove the Pregnancy Test device and pipette from the pouch and use immediately.
3. Lay the test device down on a flat, non-absorbent, clean surface.
4. Using the sample dropper, collect urine from the sample container and slowly dispense 3 drops (120µl) of urine into the cassette test well.
5. Read the result at 5 minutes. Do **not** read the result after 15 minutes.

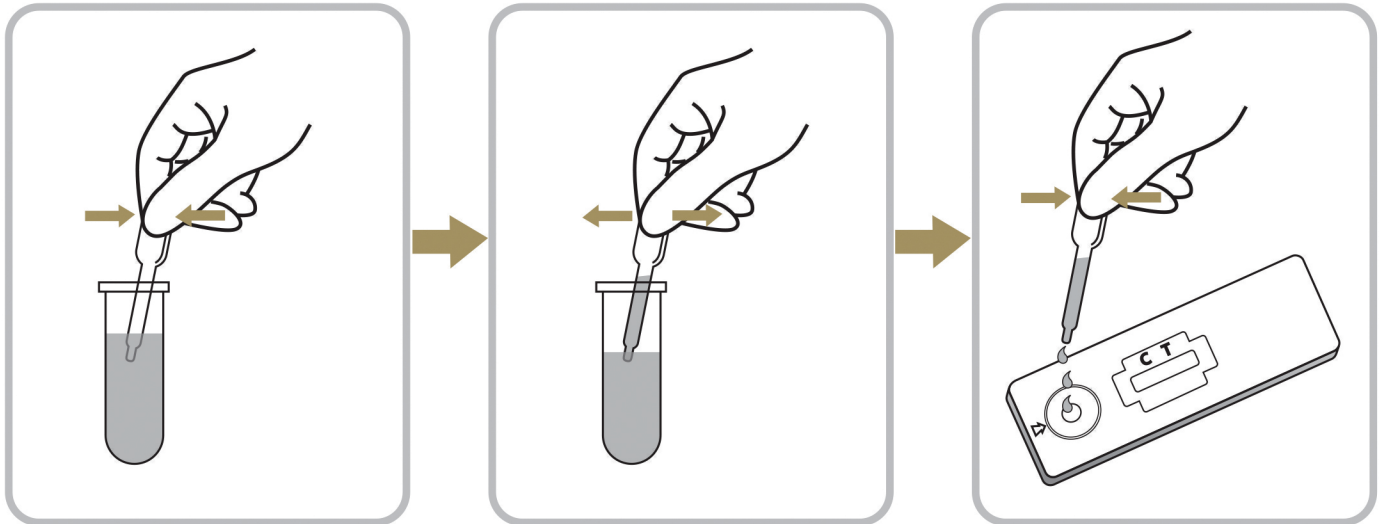